

REGAL SERIES X" x X" POWER COLUMN

FULL LENGTH (M/N: W-PCXXFL-ICU)



Install
Instructions

FOR PRESENTATION
PURPOSES ONLY
Sep-25-2025
NOT FOR
CONSTRUCTION USE

DRAWING # 00

TYPE:	QUANTITY:
PATIENT ROOM	0

IF MIRRORED UNITS ARE REQUIRED, THE
CUSTOMER MUST INDICATE THE QUANTITY
BELOW


QTY MIRRORED: _____

SYSTEM DETAILS

SYM	QTY	DESCRIPTION
PC	1	POWER COLUMN, X" x X"
FM	1	FLOOR MOUNTED SS BASE
VC	4	VERTICAL EQUIPMENT TRACK
O	1	AMICO GAS OUTLET, OXYGEN
A	1	AMICO GAS OUTLET, MED AIR
V	1	AMICO GAS OUTLET, VACUUM
E2	4	RECEPTACLE - DUPLEX RED
N2	4	RECEPTACLE - DUPLEX IVORY
NC	1	PROVISION - NURSE CALL 3G
DA	1	PROVISION - DATA 1G
MC	1	MONITOR CHANNEL - 20.5"
SLM	1	VACUUM SLIDE

IMPORTANT NOTES:

- STANDARD AVAILABLE SIZES: 9"x9", 9"x12", 12"x12", 9"x19", 19"x19", 9"x28"
- MAINTAIN A MINIMUM SPACING OF 3" C/C BETWEEN LOW VOLTAGE DEVICES, POWER OUTLETS, AND SWITCHES.
- MAINTAIN A MINIMUM SPACING OF 5" C/C BETWEEN GAS OUTLETS.
- NURSE CALL DEVICES MUST MAINTAIN A MINIMUM SPACING OF 6" C/C FROM OTHER SERVICES DUE TO THE SIZE OF THE DEVICE FACE

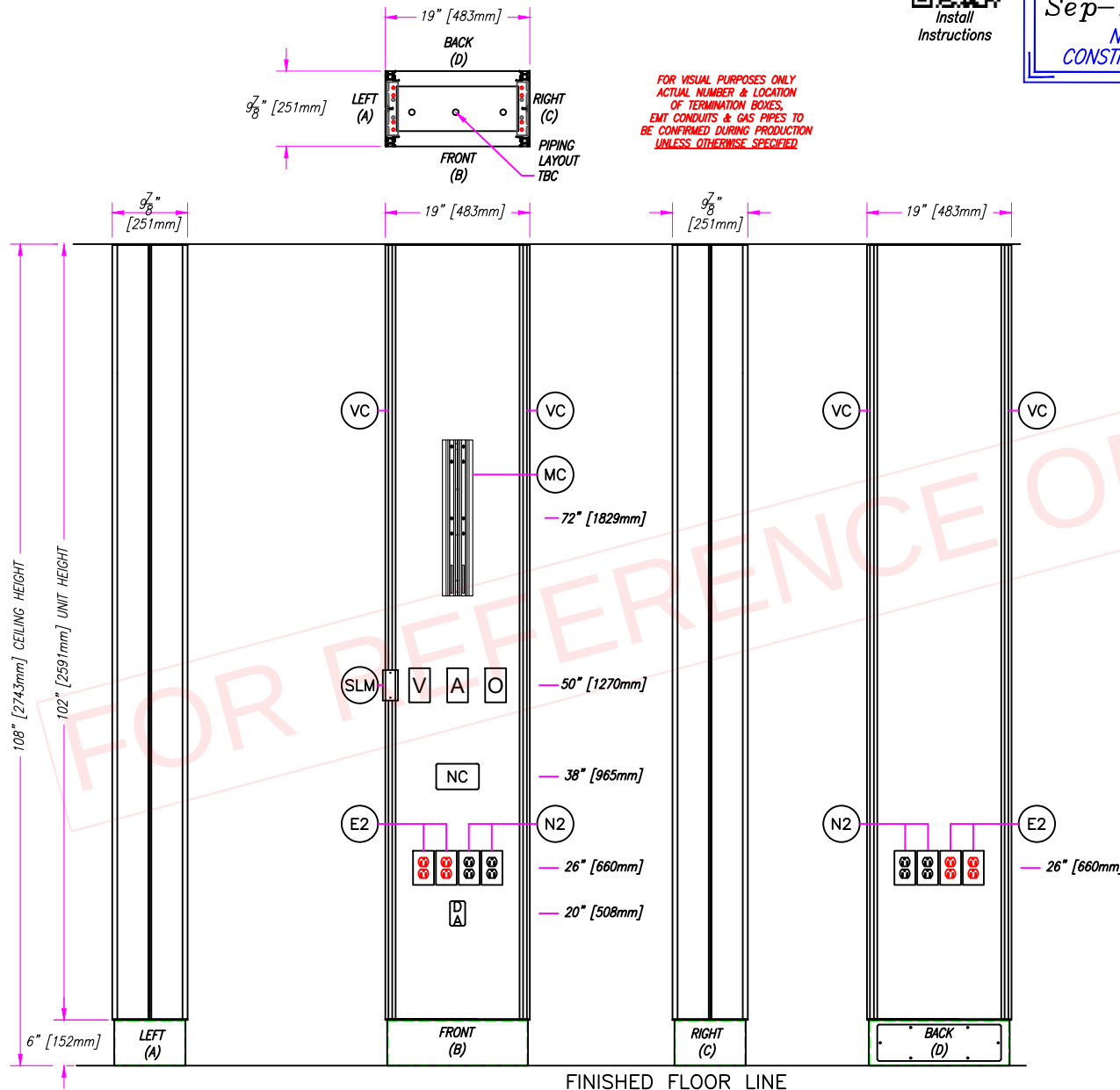

IMPORTANT: PLEASE VERIFY THAT THE ABOVE INFORMATION IS
CORRECT, AND PROVIDE THE REMAINING DETAILS.

APPROVAL - PRINT AND SIGNATURE

DATE

PHONE NO.

HOSPITAL	NURSE CALL	MFGR. : TBC	MODEL No. : TBC
LOCATION	GAS OUTLETS	MFGR. : AMICO	TYPE CONNECTION : TBC
UNIT AS SHOWN - N.T.S.	FINISH	TBC	CEILING HEIGHT TBC



85 Fulton Way
Richmond Hill, Ontario
L4B 2N4, CANADA
Toll-Free: 1-877-462-6436(T)
1-866-440-4986(F)
Tel: (905) 764-0800
Fax: (905) 764-0862

QUOTE No.
REV. BY: RJUSAYAN
PROJECT PM:
DATE: Sep-25-2025
REV. NO:


REGAL SERIES X" x X" POWER COLUMN

FULL LENGTH (M/N: W-PCXXFL-ICU)





FOR PRESENTATION
PURPOSES ONLY
Sep-25-2025
NOT FOR
CONSTRUCTION USE

$\frac{7}{8}$ "
[251mm]

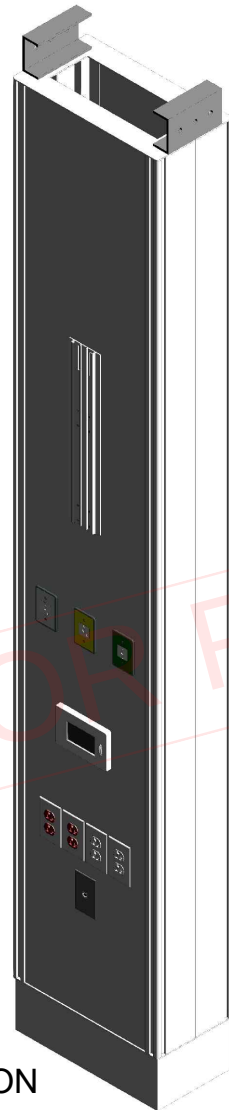

FRONT ELEVATION

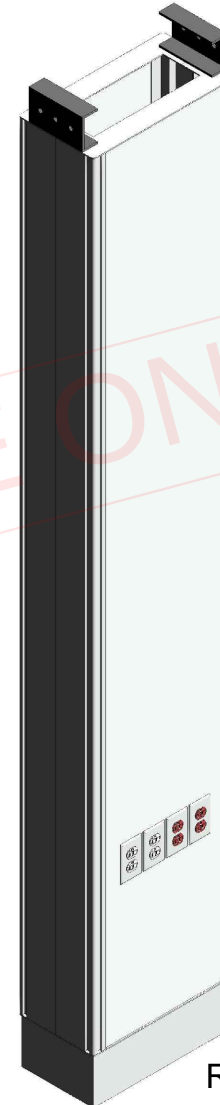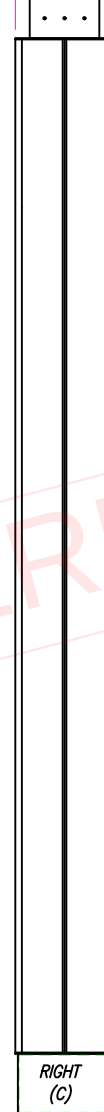

REAR ELEVATION

IMPORTANT: PLEASE VERIFY THAT THE ABOVE INFORMATION IS
CORRECT, AND PROVIDE THE REMAINING DETAILS.

APPROVAL - PRINT AND SIGNATURE

DATE

PHONE NO.

HOSPITAL	HOSPITAL	NURSE CALL	MFGR. : TBC	MODEL No. : TBC	 85 Fulton Way Richmond Hill, Ontario L4B 2N4, CANADA Toll-Free: 1-877-462-6436(T) 1-866-440-4986(F) Tel: (905) 764-0800 Fax: (905) 764-0862 www.amico.com	QUOTE No.
LOCATION	LOCATION	GAS OUTLETS	MFGR. : AMICO	TYPE CONNECTION : TBC		REV. BY: RJUSAYAN PROJECT PM: DATE: Sep-25-2025
UNIT AS SHOWN - N.T.S.		FINISH	TBC	CEILING HEIGHT		REV. NO: