

PLEASE SPECIFY (IF NOT SPECIFIED AMICO  
WILL PROVIDE DEFAULTS LISTED BELOW):

STUD HEIGHT ( MAX. 12'-0" ): TBC

STUD SIZE ( 3 5/8" TYP. ): TBC

STUD GAUGE ( #16 TYP. ): TBC

WALL THICKNESS ( 5/8" TYP. ): TBC

# REGAL SERIES RECESSED PANELWALL EXTENDED LENGTH GENERAL CARE (M/N: FWXXXX-SR-GEN-P-ZZ)



Install  
Instructions

FOR PRESENTATION  
PURPOSES ONLY  
Sep-17-2025  
NOT FOR  
CONSTRUCTION USE

## DRAWING # 00

TYPE:

PATIENT ROOM

QUANTITY:

00

IF MIRRORED UNITS ARE REQUIRED, THE  
CUSTOMER MUST INDICATE THE QUANTITY  
BELOW


QTY MIRRORED: \_\_\_\_\_

### SYSTEM DETAILS

| SYM | QTY | DESCRIPTION               |
|-----|-----|---------------------------|
| V   | 1   | AMICO GAS OUTLET, VACUUM  |
| A   | 1   | AMICO GAS OUTLET, VACUUM  |
| O   | 2   | AMICO GAS OUTLET, OXYGEN  |
| E2  | 4   | RECEPTACLE - DUPLEX RED   |
| NC  | 1   | PROVISION - NURSE CALL 3G |
| DA  | 2   | PROVISION - DATA 1G       |
| SLM | 1   | VACUUM SLIDE              |
|     | 2   | STUD ASSEMBLY             |

### IMPORTANT NOTES:

- STANDARD AVAILABLE WIDTHS: 14", 19" AND 24".
- MAINTAIN A MINIMUM SPACING OF 3" C/C BETWEEN LOW VOLTAGE DEVICES, POWER OUTLETS, AND SWITCHES.
- MAINTAIN A MINIMUM SPACING OF 5" C/C BETWEEN GAS OUTLETS.
- NURSE CALL DEVICES MUST MAINTAIN A MINIMUM SPACING OF 6" C/C FROM OTHER SERVICES DUE TO THE SIZE OF THE DEVICE FACE
- STUD SIZES AVAILABLE ON 2 1/2", 3 5/8", 4", 6" AND 8".
- BACK-TO-BACK UNITS REQUIRES MINIMUM OF 3 5/8" STUDS.
- SLOPED TOP CAP IS REQUIRED WHEN UNIT DOES NOT EXTEND TO THE CEILING.
- FLAT TOP CAP TO BE USED WHEN UNIT EXTENDS TO THE CEILING.
- DIALYSIS BOXES ARE NOT COMPATIBLE WITH THIS HEADWALL PRODUCT.

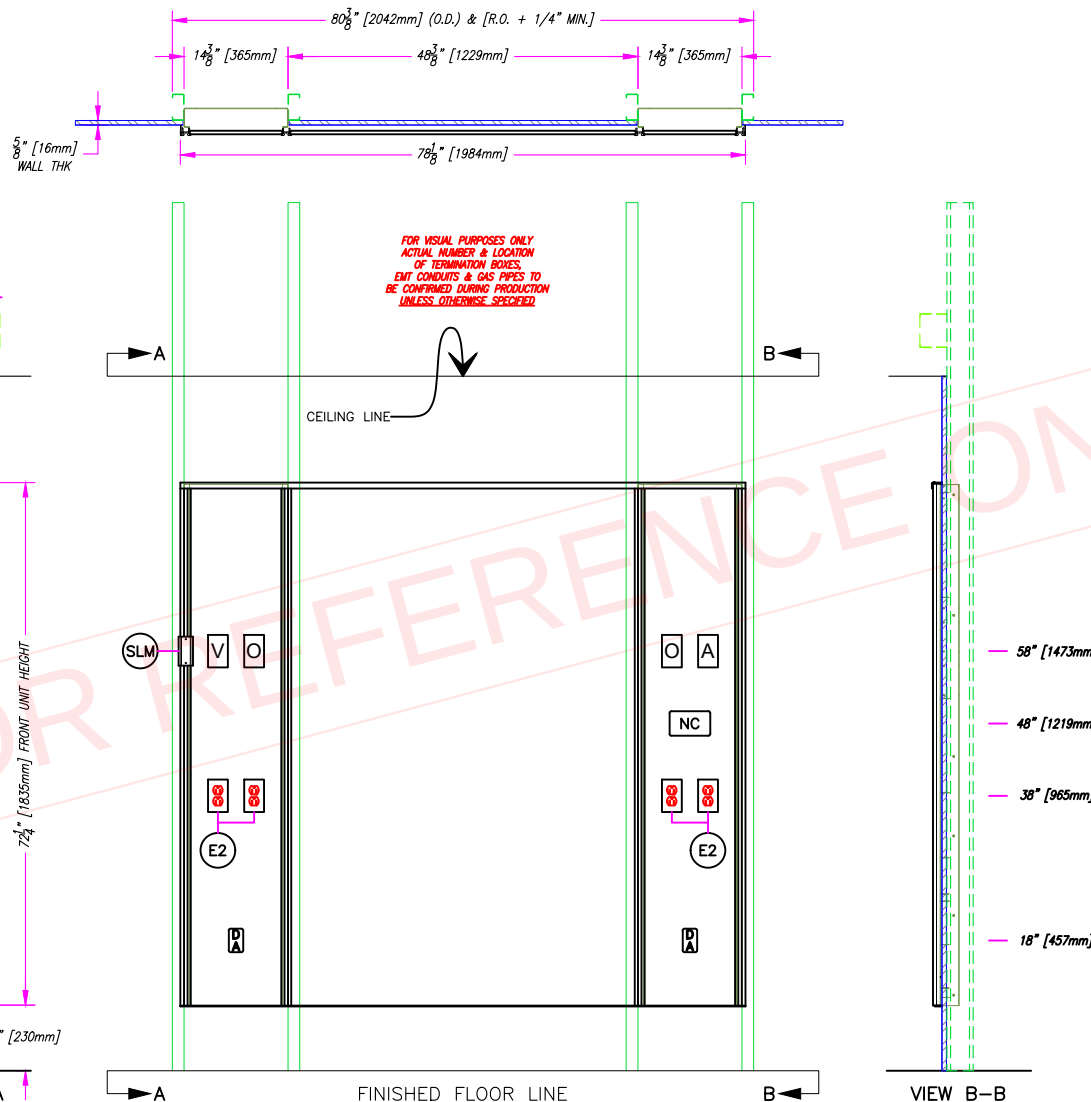

IMPORTANT: PLEASE VERIFY THAT THE ABOVE INFORMATION IS  
CORRECT, AND PROVIDE THE REMAINING DETAILS.

APPROVAL - PRINT AND SIGNATURE

DATE

PHONE NO.

|                        |          |             |               |                       |                                                                                                        |                                                                                                                                                               |                                                                                |
|------------------------|----------|-------------|---------------|-----------------------|--------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------|
| HOSPITAL               | HOSPITAL | NURSE CALL  | MFGR. : TBC   | MODEL No. : TBC       | <br>www.amico.com | 85 Fulton Way<br>Richmond Hill, Ontario<br>L4B 2N4, CANADA<br>Toll-Free: 1-877-462-6436(T)<br>1-866-440-4986(F)<br>Tel: (905) 764-0800<br>Fax: (905) 764-0862 | QUOTE No.<br>REV. BY: RJUSAYAN<br>PROJECT PM:<br>DATE: Sep-17-2025<br>REV. NO: |
| LOCATION               | LOCATION | GAS OUTLETS | MFGR. : AMICO | TYPE CONNECTION : TBC |                                                                                                        |                                                                                                                                                               |                                                                                |
| UNIT AS SHOWN - N.T.S. |          | FINISH      | TBC           | CEILING HEIGHT        |                                                                                                        |                                                                                                                                                               |                                                                                |

PLEASE SPECIFY (IF NOT SPECIFIED AMICO  
WILL PROVIDE DEFAULTS LISTED BELOW):

STUD HEIGHT ( MAX. 12'-0" ) : TBC

STUD SIZE ( 3 5/8" TYP. ) : TBC

STUD GAUGE ( #16 TYP. ) : TBC

WALL THICKNESS ( 5/8" TYP. ) : TBC

# REGAL SERIES RECESSED PANELWALL EXTENDED LENGTH GENERAL CARE (M/N: FWXXXX-SR-GEN-P-ZZ)



Install  
Instructions

FOR PRESENTATION  
PURPOSES ONLY  
Sep-17-2025  
NOT FOR  
CONSTRUCTION USE

DRAWING # 00

TYPE:

PATIENT ROOM

QUANTITY:

00



BACK TO BACK

## SYSTEM DETAILS

| SYM | QTY | DESCRIPTION               |
|-----|-----|---------------------------|
| V   | 2   | AMICO GAS OUTLET, VACUUM  |
| A   | 2   | AMICO GAS OUTLET, VACUUM  |
| O   | 4   | AMICO GAS OUTLET, OXYGEN  |
| E2  | 8   | RECEPTACLE - DUPLEX RED   |
| NC  | 2   | PROVISION - NURSE CALL 3G |
| DA  | 2   | PROVISION - DATA 1G       |
| SLM | 2   | VACUUM SLIDE              |
|     | 2   | STUD ASSEMBLY             |

### IMPORTANT NOTES:

- STANDARD AVAILABLE WIDTHS: 14", 19" AND 24".
- MAINTAIN A MINIMUM SPACING OF 3" C/C BETWEEN LOW VOLTAGE DEVICES, POWER OUTLETS, AND SWITCHES.
- MAINTAIN A MINIMUM SPACING OF 5" C/C BETWEEN GAS OUTLETS.
- NURSE CALL DEVICES MUST MAINTAIN A MINIMUM SPACING OF 6" C/C FROM OTHER SERVICES DUE TO THE SIZE OF THE DEVICE FACE
- STUD SIZES AVAILABLE ON 2 1/2", 3 5/8", 4", 6" AND 8".
- BACK-TO-BACK UNITS REQUIRES MINIMUM OF 3 5/8" STUDS.
- SLOPED TOP CAP IS REQUIRED WHEN UNIT DOES NOT EXTEND TO THE CEILING.
- FLAT TOP CAP TO BE USED WHEN UNIT EXTENDS TO THE CEILING.
- DIALYSIS BOXES ARE NOT COMPATIBLE WITH THIS HEADWALL PRODUCT.

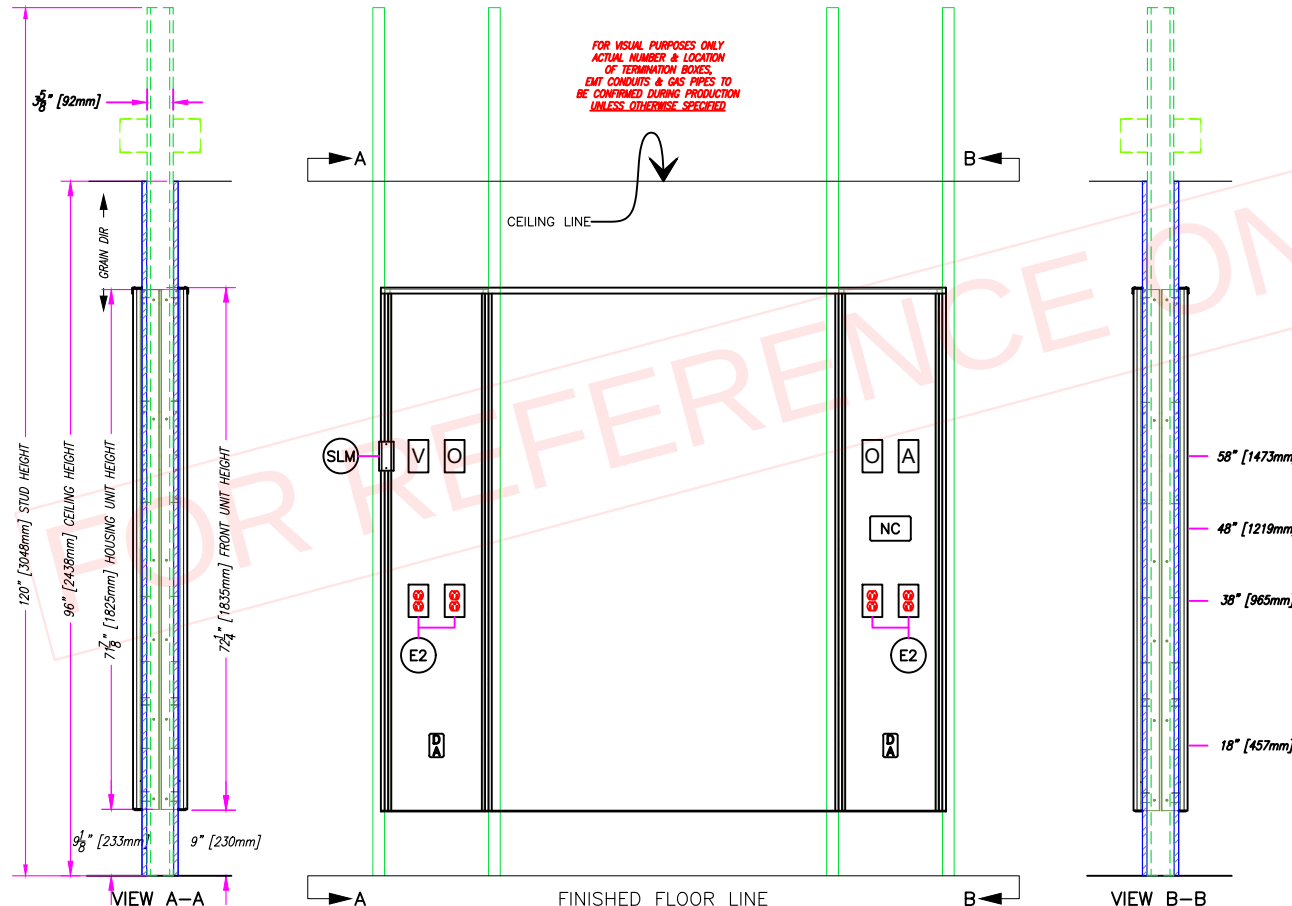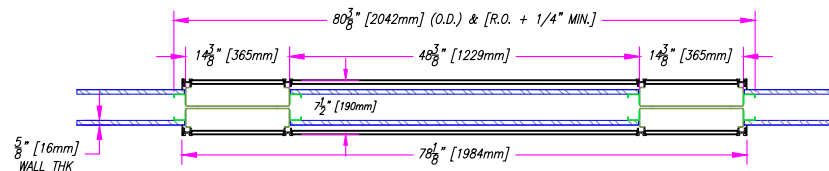

IMPORTANT: PLEASE VERIFY THAT THE ABOVE INFORMATION IS  
CORRECT, AND PROVIDE THE REMAINING DETAILS.

APPROVAL - PRINT AND SIGNATURE

DATE

PHONE NO.

|                        |          |             |               |                       |                                                                                                        |                                                                                                                                                               |                                                                                |
|------------------------|----------|-------------|---------------|-----------------------|--------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------|--------------------------------------------------------------------------------|
| HOSPITAL               | HOSPITAL | NURSE CALL  | MFGR. : TBC   | MODEL No. : TBC       | <br>www.amico.com | 85 Fulton Way<br>Richmond Hill, Ontario<br>L4B 2N4, CANADA<br>Toll-Free: 1-877-462-6436(T)<br>1-866-440-4986(F)<br>Tel: (905) 764-0800<br>Fax: (905) 764-0862 | QUOTE No.<br>REV. BY: RJUSAYAN<br>PROJECT PM:<br>DATE: Sep-17-2025<br>REV. NO: |
| LOCATION               | LOCATION | GAS OUTLETS | MFGR. : AMICO | TYPE CONNECTION : TBC |                                                                                                        |                                                                                                                                                               |                                                                                |
| UNIT AS SHOWN - N.T.S. |          | FINISH      | TBC           | CEILING HEIGHT        |                                                                                                        |                                                                                                                                                               |                                                                                |

PLEASE SPECIFY (IF NOT SPECIFIED AMICO  
WILL PROVIDE DEFAULTS LISTED BELOW):

STUD HEIGHT ( MAX. 12'-0" ): TBC

STUD SIZE ( 3 5/8" TYP. ): TBC

STUD GAUGE ( #16 TYP. ): TBC

WALL THICKNESS ( 5/8" TYP. ): TBC

# REGAL SERIES RECESSED PANELWALL EXTENDED LENGTH GENERAL CARE (M/N: FWXXXX-SR-GEN-P-ZZ)




Install  
Instructions

FOR PRESENTATION  
PURPOSES ONLY  
Sep-17-2025  
NOT FOR  
CONSTRUCTION USE

## BACK TO BACK UNIT DETAIL



B2B AS SHOWN SIDE

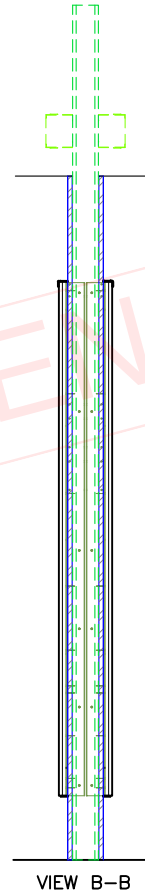

B2B MIRRORED SIDE

IMPORTANT: PLEASE VERIFY THAT THE ABOVE INFORMATION IS  
CORRECT, AND PROVIDE THE REMAINING DETAILS.

APPROVAL - PRINT AND SIGNATURE

DATE

PHONE NO.

|                        |          |             |               |                       |                                                                                                        |                                                                                                                                                               |                                                                                    |
|------------------------|----------|-------------|---------------|-----------------------|--------------------------------------------------------------------------------------------------------|---------------------------------------------------------------------------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------|
| HOSPITAL               | HOSPITAL | NURSE CALL  | MFGR. : TBC   | MODEL No. : TBC       | <br>www.amico.com | 85 Fulton Way<br>Richmond Hill, Ontario<br>L4B 2N4, CANADA<br>Toll-Free: 1-877-462-6436(T)<br>1-866-440-4986(F)<br>Tel: (905) 764-0800<br>Fax: (905) 764-0862 | QUOTE No.<br><br>REV. BY: RJUSAYAN<br>PROJECT PM:<br>DATE: Sep-17-2025<br>REV. NO: |
| LOCATION               | LOCATION | GAS OUTLETS | MFGR. : AMICO | TYPE CONNECTION : TBC |                                                                                                        |                                                                                                                                                               |                                                                                    |
| UNIT AS SHOWN - N.T.S. |          | FINISH      | TBC           | CEILING HEIGHT        |                                                                                                        |                                                                                                                                                               |                                                                                    |
|                        |          |             |               | TBC                   |                                                                                                        |                                                                                                                                                               |                                                                                    |