

PLEASE SPECIFY (IF NOT SPECIFIED AMICO
WILL PROVIDE DEFAULTS LISTED BELOW):

STUD HEIGHT (MAX. 12'-0"): TBC

STUD SIZE (3 5/8" TYP.): TBC

STUD GAUGE (#16 TYP.): TBC

WALL THICKNESS (5/8" TYP.): TBC

REGAL SERIES RECESSED FLATWALL

[EXTENDED] LENGTH GENERAL CARE

[W/ STUD SINGLE SIDE]

(M/N: FWXXXX-RC-GEN-P)



Install
Instructions

FOR PRESENTATION
PURPOSES ONLY
Sep-17-2025
NOT FOR
CONSTRUCTION USE

DRAWING # 00

TYPE:

PATIENT ROOM

QUANTITY:

00

IF MIRRORED UNITS ARE REQUIRED, THE
CUSTOMER MUST INDICATE THE
QUANTITY BELOW


QTY MIRRORED: _____

SYSTEM DETAILS


SYM	QTY	DESCRIPTION
V	1	AMICO GAS OUTLET, VACUUM
O	1	AMICO GAS OUTLET, OXYGEN
E2	2	RECEPTACLE - DUPLEX RED
NC	1	PROVISION - NURSE CALL 3G
DA	1	PROVISION - DATA 1G
SLM	1	VACUUM SLIDE
	1	STUD ASSEMBLY

IMPORTANT NOTES:

- STANDARD AVAILABLE WIDTHS: 14", 19" AND 24".
- MAINTAIN A MINIMUM SPACING OF 3" C/C BETWEEN LOW VOLTAGE DEVICES, POWER OUTLETS, AND SWITCHES.
- MAINTAIN A MINIMUM SPACING OF 5" C/C BETWEEN GAS OUTLETS.
- NURSE CALL DEVICES MUST MAINTAIN A MINIMUM SPACING OF 6" C/C FROM OTHER SERVICES DUE TO THE SIZE OF THE DEVICE FACE
- STUD SIZES AVAILABLE ON 2 1/2", 3 5/8", 4", 6" AND 8".
- BACK-TO-BACK UNITS REQUIRES MINIMUM OF 3 5/8" STUDS.
- SLOPED TOP CAP IS REQUIRED WHEN UNIT DOES NOT EXTEND TO THE CEILING.
- FLAT TOP CAP TO BE USED WHEN UNIT EXTENDS TO THE CEILING.
- DIALYSIS BOXES ARE NOT COMPATIBLE WITH THIS HEADWALL PRODUCT.

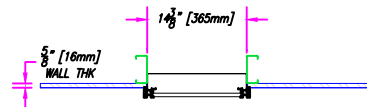

FOR VISUAL PURPOSES ONLY
ACTUAL NUMBER & LOCATION
OF TERMINATION BOXES,
EMT CONDUITS & GAS PIPES TO
BE CONFIRMED DURING PRODUCTION
UNLESS OTHERWISE SPECIFIED

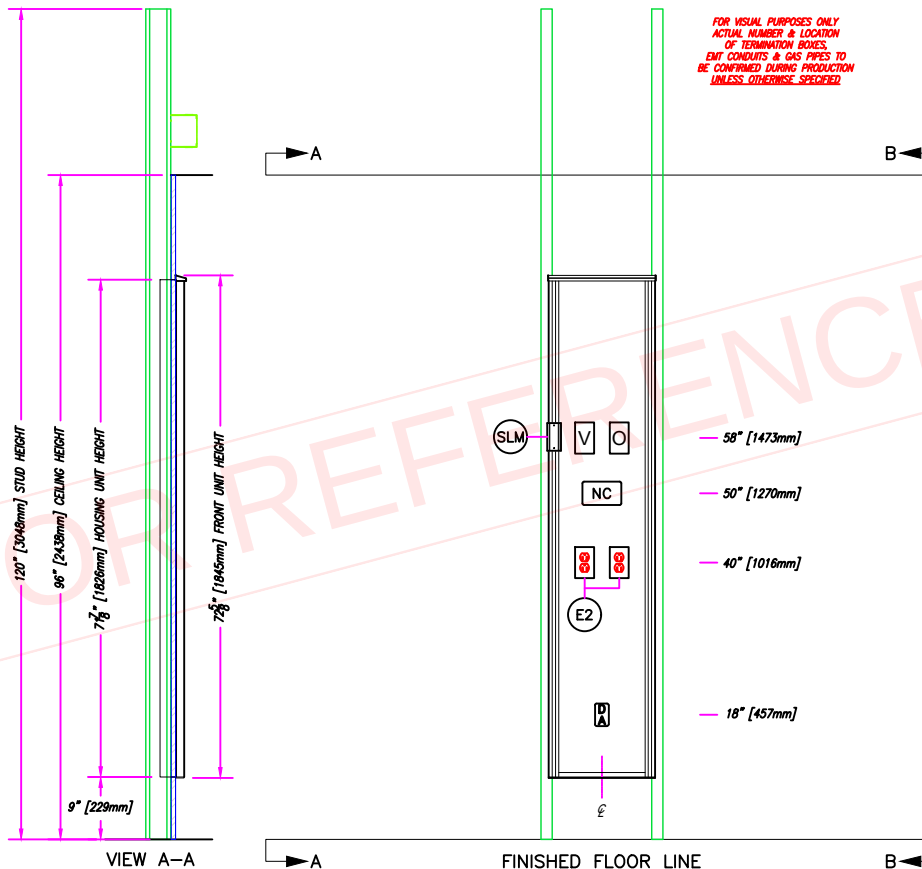

IMPORTANT: PLEASE VERIFY THAT THE ABOVE INFORMATION IS
CORRECT, AND PROVIDE THE REMAINING DETAILS.

APPROVAL - PRINT AND SIGNATURE

DATE

PHONE NO.

HOSPITAL	HOSPITAL	NURSE CALL	MFGR. : TBC	MODEL No. : TBC
LOCATION	LOCATION	GAS OUTLETS	MFGR. : AMICO	TYPE CONNECTION : TBC
	UNIT AS SHOWN - N.T.S.	FINISH	TBC	CEILING HEIGHT TBC



85 Fulton Way
Richmond Hill, Ontario
L4B 2N4, CANADA
Toll-Free: 1-877-462-6436(T)
1-866-440-4986(F)
Tel: (905) 764-0800
Fax: (905) 764-0862

QUOTE No.

REV. BY: RJUSAYAN
PROJECT PM:
DATE: Sep-17-2025
REV. NO:

PLEASE SPECIFY (IF NOT SPECIFIED AMICO
WILL PROVIDE DEFAULTS LISTED BELOW):

STUD HEIGHT (MAX. 12'-0") : TBC

STUD SIZE (3 5/8" TYP.) : TBC

STUD GAUGE (#16 TYP.) : TBC

WALL THICKNESS (5/8" TYP.) : TBC

REGAL SERIES RECESSED FLATWALL

[EXTENDED] LENGTH GENERAL CARE

[W/ STUD BACK TO BACK],

(M/N: FWXXXX-BB-GEN-P)



Install
Instructions

FOR PRESENTATION
PURPOSES ONLY
Sep-17-2025
NOT FOR
CONSTRUCTION USE

DRAWING # 00

TYPE:

PATIENT ROOM

QUANTITY:

00



BACK TO BACK

SYSTEM DETAILS

SYM	QTY	DESCRIPTION
V	2	AMICO GAS OUTLET, VACUUM
O	2	AMICO GAS OUTLET, OXYGEN
E2	4	RECEPTACLE - DUPLEX RED
NC	2	PROVISION - NURSE CALL 3G
DA	2	PROVISION - DATA 1G
SLM	2	VACUUM SLIDE
	1	STUD ASSEMBLY

IMPORTANT NOTES:

- STANDARD AVAILABLE WIDTHS: 14", 19" AND 24".
- MAINTAIN A MINIMUM SPACING OF 3" C/C BETWEEN LOW VOLTAGE DEVICES, POWER OUTLETS, AND SWITCHES.
- MAINTAIN A MINIMUM SPACING OF 5" C/C BETWEEN GAS OUTLETS.
- NURSE CALL DEVICES MUST MAINTAIN A MINIMUM SPACING OF 6" C/C FROM OTHER SERVICES DUE TO THE SIZE OF THE DEVICE FACE
- STUD SIZES AVAILABLE ON 2 1/2", 3 5/8", 4", 6" AND 8".
- BACK-TO-BACK UNITS REQUIRES MINIMUM OF 3 5/8" STUDS.
- SLOPED TOP CAP IS REQUIRED WHEN UNIT DOES NOT EXTEND TO THE CEILING.
- FLAT TOP CAP TO BE USED WHEN UNIT EXTENDS TO THE CEILING.
- DIALYSIS BOXES ARE NOT COMPATIBLE WITH THIS HEADWALL PRODUCT.

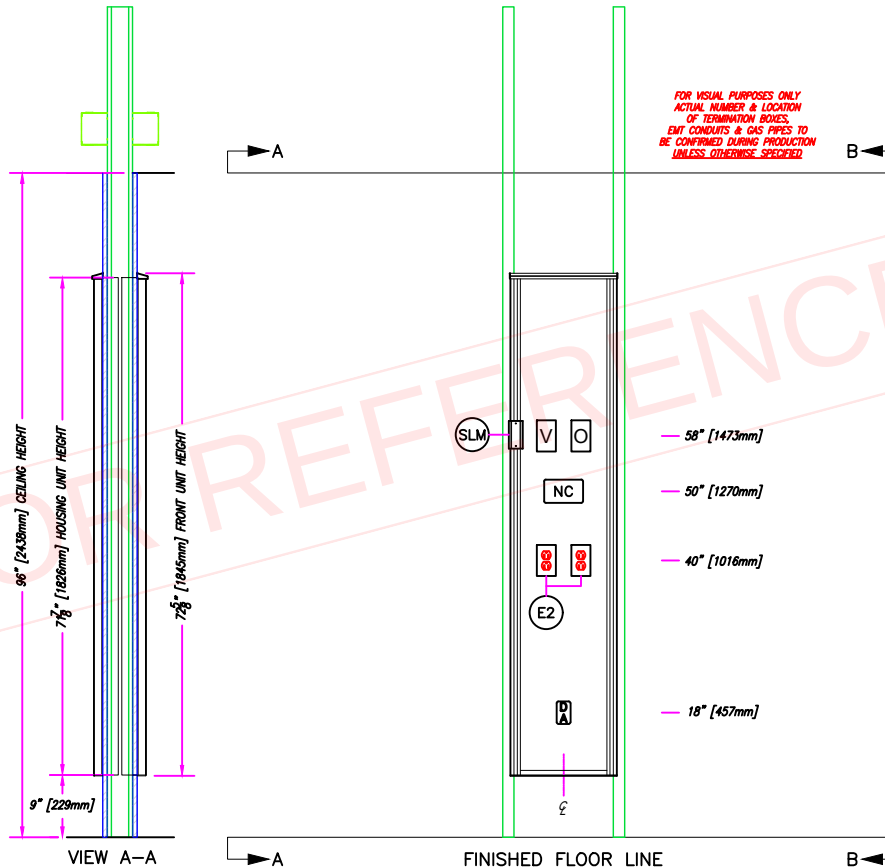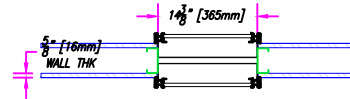

FOR VISUAL PURPOSES ONLY
ACTUAL NUMBER & LOCATION
OF TERMINATION BOXES,
EXIT CONDUITS & GAS PIPES TO
BE CONFIRMED DURING PRODUCTION
UNLESS OTHERWISE SPECIFIED

IMPORTANT: PLEASE VERIFY THAT THE ABOVE INFORMATION IS
CORRECT, AND PROVIDE THE REMAINING DETAILS.

APPROVAL - PRINT AND SIGNATURE

DATE

PHONE NO.

HOSPITAL	HOSPITAL	NURSE CALL	MFGR. : TBC	MODEL No. : TBC
LOCATION	LOCATION	GAS OUTLETS	MFGR. : AMICO	TYPE CONNECTION : TBC
	UNIT AS SHOWN - N.T.S.	FINISH	TBC	CEILING HEIGHT TBC



85 Fulton Way
Richmond Hill, Ontario
L4B 2N4, CANADA
Toll-Free: 1-877-462-6436(T)
1-866-440-4986(F)
Tel: (905) 764-0800
Fax: (905) 764-0862

QUOTE No.

REV. BY: RJUSAYAN
PROJECT PM:
DATE: Sep-17-2025
REV. NO:

PLEASE SPECIFY (IF NOT SPECIFIED AMICO
WILL PROVIDE DEFAULTS LISTED BELOW):

STUD HEIGHT (MAX. 12'-0"): TBC

STUD SIZE (3 5/8" TYP.): TBC

STUD GAUGE (#16 TYP.): TBC

WALL THICKNESS (5/8" TYP.): TBC

REGAL SERIES RECESSED FLATWALL

[EXTENDED] LENGTH GENERAL CARE

[W/ STUD BACK TO BACK]


(M/N: FWXXXX-BB-GEN-P)




Install
Instructions

FOR PRESENTATION
PURPOSES ONLY
Sep-17-2025
NOT FOR
CONSTRUCTION USE


BACK TO BACK
UNIT DETAIL

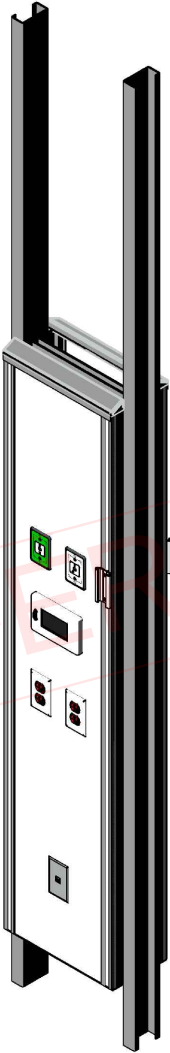

B2B AS SHOWN SIDE

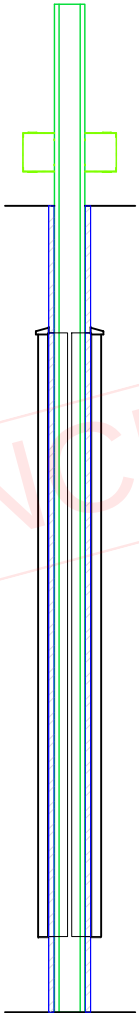

VIEW A-A

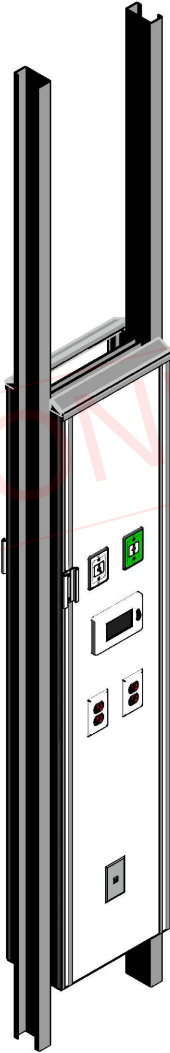

B2B MIRRORED SIDE

IMPORTANT: PLEASE VERIFY THAT THE ABOVE INFORMATION IS
CORRECT, AND PROVIDE THE REMAINING DETAILS.

APPROVAL - PRINT AND SIGNATURE

DATE

PHONE NO.

HOSPITAL	HOSPITAL
LOCATION	LOCATION
UNIT AS SHOWN - N.T.S.	

NURSE CALL	MFGR. : TBC	MODEL No. : TBC
GAS OUTLETS	MFGR. : AMICO	TYPE CONNECTION : TBC
FINISH	TBC	CEILING HEIGHT TBC



85 Fulton Way
Richmond Hill, Ontario
L4B 2N4, CANADA
Toll-Free: 1-877-462-6436(T)
1-866-440-4986(F)
Tel: (905) 764-0800
Fax: (905) 764-0862

QUOTE No.
REV. BY: RJUSAYAN
PROJECT PM:
DATE: Sep-17-2025
REV. NO: