

Anesthesia Mounting Solutions



Amico

Anesthesia Mounting Products

Adjustable Height AHM Monitor Arm

Features

- 180° side to side motion
 - Ability to lock arm in any position
 - Can be mounted directly to any equipment rail
 - Integrated cable management
 - 17" (43.2 cm) of height adjustment
 - 30 lbs (14 kg) arm rating - heavy duty option available for loads up to 40 lbs (18.2 kg)
 - Internal gas spring provides counterbalance
 - Durable, aluminum body
 - Quick release pin for single hand monitor removal
-



Side to Side SSM Monitor Arm

Features

- 180° side to side motion
 - Can be mounted directly to any equipment rail
 - Integrated cable management
 - 9" or 12" (23 cm or 30 cm) single and dual arm lengths available
 - Tilt and swivel adjustments
 - 30 lbs (13.6 kg) single arm rating
 - 15 lbs (6.8 kg) extension arm rating
-



Adjustable Height AHM Monitor Arm

Features

- 6.5" (16.51 cm) Independent LCD height adjustment
- Entire system easily folds back less than 9" (22.9 cm) for storage when not in use
- Includes both 75/100 VESA mounting options with easy installation slots
- Integrated cable management
- Negative tilt keyboard adjustment
- Easy access to barcode scanner and mouse holders
- 17" (43.18 Cm) Overall vertical adjustment range accommodates users of any height
- Locking feature when height has been adjusted
- Counterbalanced to provide effortless movement



Anesthesia Mounting Features



Cable Management

Discrete and sleek cable management maintains a tidy workspace and encourages excellent infection control.



Custom Solutions

A variety of options are available to integrate multiple monitors, an IT charting station or additional peripherals. We work with your facility to accommodate specific needs and unique requirements.




Pre-Assembled

All of our products come pre-assembled to save you the time and cost that it takes to assemble the units. Simply slide the dovetail or MRS adapter into the track on the machine.



Improve Flexibility

Increase the range of motion for your IT charting station or patient monitor by adding an extension, but don't compromise the stability. Extensions come in a variety of sizes from 6" to 19" (15.24 cm to 48.26 cm).

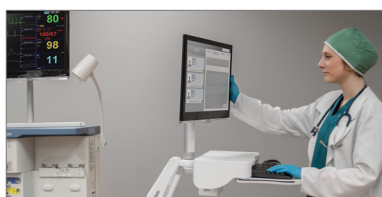

Easy to Use

Adjust the IT charting station or patient monitor with one hand while working! All of our height adjustable units are counterbalanced to improve workflow and allow the clinician to focus on patient care.

Monitor Compatibility

- Advantech
- Alaris
- Analogic
- CAS Medical
- Capsule
- Cerner
- Colin
- Criticare
- Datascope
- Datex-Ohmeda
- Dräger Medical
- Elo Touch
- Fukuda Denshi
- GE Healthcare
- Hoana Medical
- INVIVO
- Ivy Biomedical Systems
- Larsen & Toubro
- Masimo
- Mediana
- Medtronic
- Mennen Medical
- Mindray
- Nellcor
- Nihon Kohden
- Nonin
- Oridion
- Spacelabs Healthcare
- Schiller
- Spacelabs Medical
- Welch Allyn
- WITT
- ZOLL

Options available for other manufacturers, please specify model when ordering

Application Examples

Add Patient Monitoring

- Patient monitor
- Recorder unit
- Connector
- Height adjustability
- Swivel adjustment
- Cable management



Add IT Stations

- IT station for charting
- CPU Mount
- Height adjustability
- Swivel adjustment
- Portrait orientation for LCD monitor
- Cable management
- Extra writing surface



Anesthesia Machine Compatibility

- | | | |
|------------------------|----------------------------|---|
| • GE Avance | • Dräger Julian | • Datex Ohmeda Avance |
| • GE Aisys | • Dräger Tiro | • Datex Ohmeda Aisys |
| • GE Aestive | • Dräger Fabius OS | • Datex Ohmeda Aestiva |
| • GE Aespire | • Dräger Fabius MRI | • Datex Ohmeda Aespire |
| • Dräger Primus | • Dräger Fabius GS | • Datex Ohmeda ADU |
| • Dräger Narkomed 6000 | • Dräger Fabius GS Premium | • Custom - talk to a product expert today! Call 1.877.264.2697 or Email accessories@amico.com |
| • Dräger Narkomed 6400 | • Dräger Cato | |
| • Dräger Narkomed GS | • Dräger Apollo | |