

Stretchers



Titan Transport

Amico

Designed for Caregivers to Provide More Care with Less Effort

The Amico Titan Stretcher uses an innovative design and incorporates distinctive features which bring efficiency to the medical and surgical environment.

The movement of the Titan Stretcher is controlled by the caregiver with pedals located around the center of the base. Pedals are located on each side of the stretcher where one pedal is used to raise the stretcher and the other pedal is used to either lower or set the stretcher in Trendelenburg or reverse Trendelenburg position.

The Titan Stretcher also features several standard functions such as a central brake locking system, adjustable fowler section, radiolucent board, adjustable knee section, and roller bumpers.



Mobility

Whether it is repositioning the stretcher within the room, or transferring the patient to a different ward, the caregiver must be able to maneuver a stretcher safely and easily. We recognize this need and have incorporated such features in the Amico Titan Stretchers.



Ergonomic Push Handles

- Simple lift-drop design for quick, one-handed activation.
- Contoured rubber grips reduce strain and pressure on the caregiver's hands and arms.



Fifth Wheel Steering

- Enhanced maneuverability when turning corners and on straightaways.



8" Wheel Casters

- Rugged 8" (20 cm) design for added durability and ease of movement over objects such as door sills.

Options

Custom configurations available for customer requirements.



4-Caster Brake System

- A braking system that provides an upgrade to the standard Head/Foot end brake pedal system. Brake pedals are located in front of each caster on the stretcher which allows for convenient braking from all sides.



4-Sided Brake System

- A braking system that provides the ultimate in caregiver convenience. Brake pedals are located at the Head, Foot, and each Side of the stretcher, providing the most efficient braking option for the caregiver.

Accessories



Foldable IV Pole

- Permanently fixed IV pole eliminates the chance of the IV pole being misplaced.
- Folds away for quick storage with only one hand.



Chair Position

- Chair position to sit the patient up and get them mobilized quicker.

Simplicity

The Amico Titan Stretcher controls are designed to be easy-to-use and promote efficient patient care.



Steering System

- Enhanced maneuverability when turning corners and on straightaways.



Hi-lo Pump Pedal

- Pump pedals on each side of the stretcher allow the caregiver to adjust the height without leaving the patient's side.



Knee Gatch

- Allows the caregiver to adjust the angle of the Knee Section for added patient comfort or treatment positions (e.g., supine).



Foldable Siderails with Bumper

- Easy-to-use design allows for one-handed lifting with an automatic locking latch.
- Full length coverage provides increased protection against patient falls.



Radiolucent Fowler Board

- Offers caregivers the convenience of performing upper torso X-rays without removing the patient from the stretcher.



Oxygen Tank Holder and Patient Item Storage

- Large storage area provides a convenient location for all of the patient's belongings.
- Compartments for two oxygen tanks.



Down/Trendelenburg/ Reverse Trendelenburg Control Pedals

- Control Pedals on each side of the stretcher provide caregivers with full adjustability of the stretcher's height and Trendelenburg angle, without leaving the patient's side.



Brake & Neutral Pedal

- Conveniently located at the Head end and Foot end for quick access during transport.
- Durable design with highly-visible, colored inserts to differentiate between Brake and Drive.



CPR/Head Section Articulation

- Gas spring-release mechanism also allows caregivers to quickly lower the Head Section in an emergency CPR situation.
- Dual gas spring-lift mechanism for assisted movement of the Head section.
- Auto-locking allows the Head Section to be positioned at any desired angle.



Amico Therapeutic Solutions

BodyZone Foam Mattresses



Standard Dimensions

- 76" x 25" x 4" (193 cm x 63.5 cm x 10.16 cm)

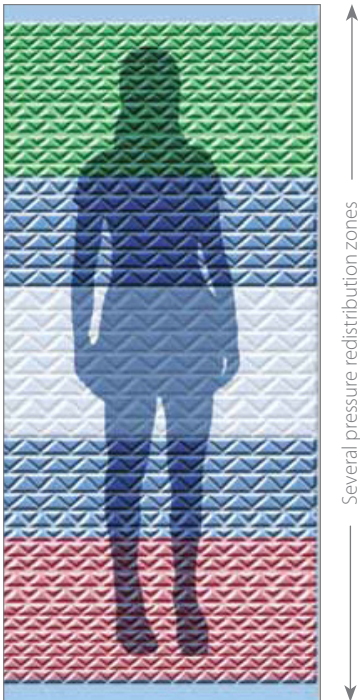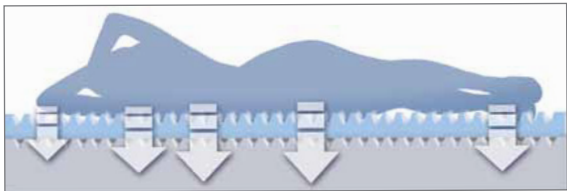

The Amico Family of Pressure Redistribution Solutions

Available for Stretchers, Beds, and Premium Beds

	50 Stretcher	100 Standard	500 Premium
VPF™	✓	✓	✓
SMT™	✓	✓	✓
Natural Ingredients	✓	✓	✓
360° Stretch Cover	✓	✓	✓
Open Cell Visco			✓

The Most Advanced Pressure Ulcer Prevention Available

- Engineered specifically for high risk areas including the head, shoulders, hips, legs, and heels.
- Envelopment zones strategically combine with body contouring zones to enhance the pressure redistribution by up to 40%.
- 360° Stretch Cover moves in every direction to minimize the forces of shear and friction.
- Patented open cell memory foam technology is proven to remove body heat at a greater rate than conventional foams and provide a cooler, drier experience.
- Our unique variable pressure foaming process provides the greatest consistency and highest quality of medical mattresses.
- Virtually zero-emissions manufacturing process.
- Patented, enclosed system that exceeds clean air act requirements.

Technical Specifications

Model Number & Approval

Model Number	S-AM1-3000
Agency Approval	cULus

Overall Dimensions

Overall Length	84.5" (214 cm)
Overall Width	33" (84 cm)

Litter Positioning

Trendelenburg/Reverse Trendelenburg	16°
Maximum Fowler Section Inclination	80° ± 5
Maximum Knee Section Inclination	35° ± 5

Height Range

Lowest	23" (58 cm)
Highest	35" (89 cm)

Surface Dimensions

Surface Width	25" (63.5 cm)
Surface Length	75" (191 cm)
Surface Thickness	4" (10 cm)

Weight Capacity & Caster Dimension

Maximum Weight Limit	550 lbs (250 kg)
Caster Diameter	8" (20 cm)

Optional Features & Accessories:

- 4"/5" (10 cm/12.7 cm) high density foam mattress with or without fire barrier
- Utility tray
- Foley bag holder
- Overbed tables
- Permanent IV pole

Amico Beds Corporation reserves the right to change product specifications without notice.

Warranty: One year parts only.

Amico Beds Corporation | 85 Fulton Way, Richmond Hill, ON L4B 2N4, Canada
Toll Free Tel: 1.877.462.6426 | Tel: 905.764.0800 | Fax: 905.764.0862
Email: info@amico.com | www.amico.com

AB-BR-TITAN 12.13.2022



1 YEAR
warranty

