

# Nova

## Motorized Exam LED Light



Amico



# The Nova



The **Nova** provides targeted exam lighting without the clutter. It is a fully recessed, motorized exam light that offers a convenient alternative to ceiling-mounted exam lighting. Engineered with 360 degrees of rotation and surgical-grade LED optics, the Nova provides high-level illumination to all corners of the patient bed. This fixture can be mounted anywhere within the ceiling to meet the needs of clinicians in all areas of the hospital environment.

Labor & Delivery

Patient Rooms

Intensive Care Units

Emergency Departments

Cardiac Cath Lab

Hybrid Operating Room







# Convenient Capacitive Touch Keypad

The **Nova Wall Control Keypad** offers precise and convenient control of light intensity and direction. Adjust up to two fixtures on one keypad with easy-to-use capacitive touch controls. Its compact design fits seamlessly into a single-gang space and can be mounted anywhere in the room. For infection control, the Nova keypad is designed with a Mylar cover and can withstand any hospital grade disinfectant.

Multiple keypads can be installed in each room.





# Clutter-Free Design

The Nova can replace any ceiling or wall-mounted exam lighting within the hospital environment. Its compact and discrete design **provides a clean patient room aesthetic.** This recessed fixture is ideal for use in low-ceiling environments and provides a **clutter-free atmosphere to enhance patient comfort.**



The minimalist design of the Nova Keypad fits seamlessly into any single-gang space. Mount the low-voltage keypad anywhere within the patient room without disrupting the visual aesthetic or taking away space for critical healthcare equipment.





# Maximize Clinical Performance

Customize your Nova configuration to maximize light output onto the patient bed. Surgical grade optics provide advanced visual acuity and superior shadow control –while 360 degrees of rotation offers clinicians increased functionality. Easily adjust light direction to provide targeted exam lighting to all areas of the patient bed — and beyond.

- **360 Degrees** of rotation provides targeted lighting to the patient bed for increased functionality
- **Compact, clutter-free design** provides a visual comfort to the patient
- **Fully recessed** exam light reduces potential hazards in low-ceiling environments and provides more space for critical healthcare equipment
- **55,000 Lux of light at task** delivers high-level illumination to support clinicians with advanced visual acuity
- **Multiple ceiling configurations** increase functionality and ensures maximum light output onto the patient bed
- **A low-voltage capacitive touch keypad** offers easy-to-use controls that ensure user safety






Features:

- 5-year limited warranty
- 15 LED optics
- 50,000 hour LED life span
- D<sub>10</sub> - 9.1°
- Motor life expectancy: 500,000 cycles
- Dimming: 1%
- IP65 Ingress protection
- Fully recessed installation

Ceiling Chart:

Ceiling Height (A)	Light Output in LUX at Nadir (B)	Light Output in LUX at Extreme Position (C)	Maximum Setback (D)
10'	55,000	21,500	5'6"
9'6"	60,700	25,000	5'1"
9'	66,300	28,800	4'8"
8'6"	72,700	32,500	4'4"
8'	80,700	36,200	3'11"

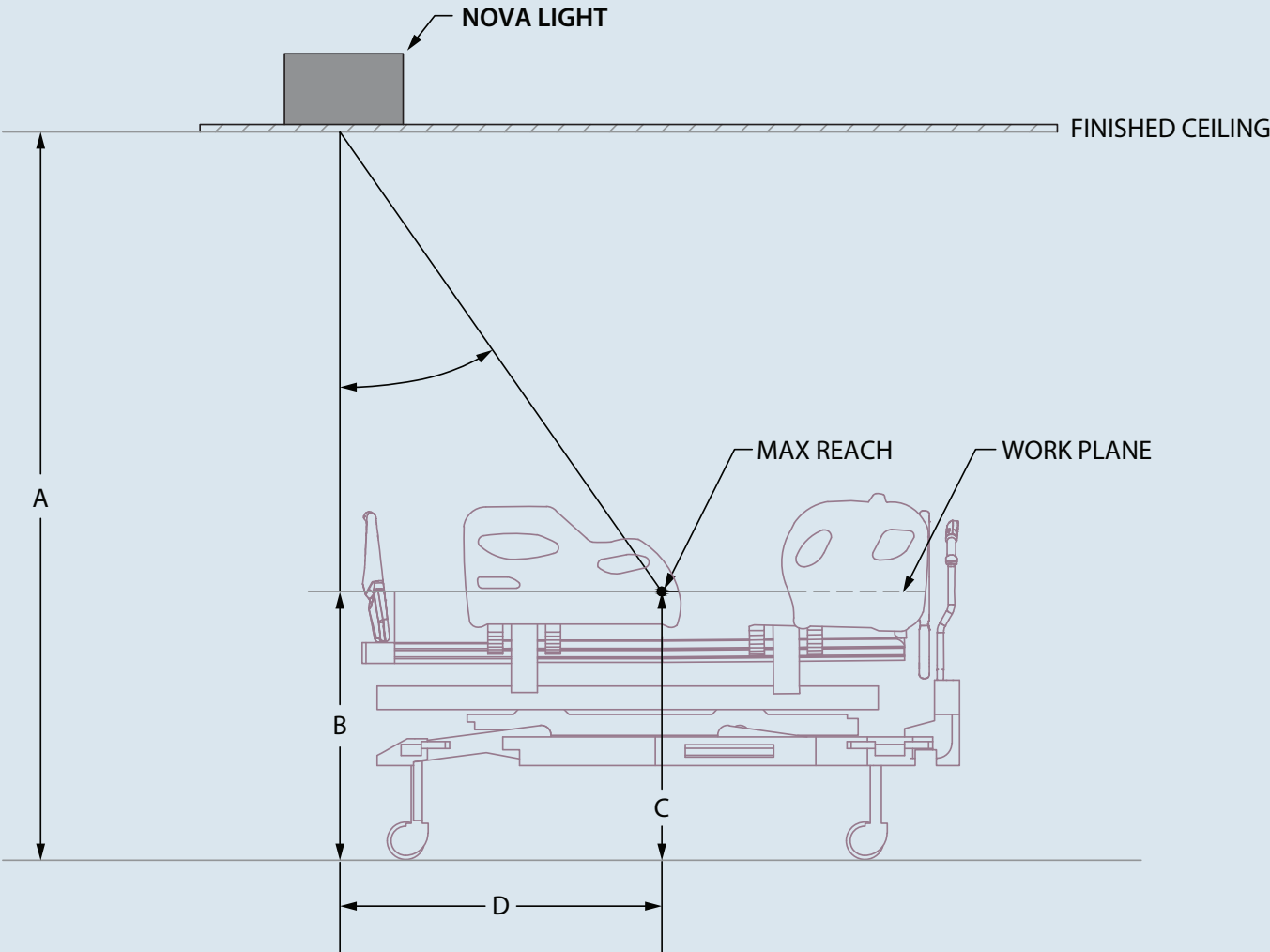

Amico Lights Corporation  
122-B East Beaver Creek Road, Richmond Hill, ON L4B 1G6, Canada  
Toll Free Phone: 1.877.462.6426 | Toll Free Fax: 1.866.440.4986  
Tel: 905.764.0800 | Fax: 905.764.0862  
Email: alt-sales@amico.com | www.amico.com

Visit **amico.com** for more information  
about our **Nova Motorized Exam Light**

