

Vega 80 LED

Exam Light



*A*mico



Power. Clarity. Efficiency.

Exam lighting you can rely on when the stakes are high.

The **Vega 80 LED Exam Light** delivers 80,000 LUX of high-intensity illumination to perform minor procedures and examinations with accuracy. The lightweight and compact design make it easy to position the fixture for maximal efficiency. So you can stay focused on the task at hand.

- 18 precisely positioned LED optics deliver high-level illumination and superior shadow control
- Color temperature of 4500K and a CRI of 95+ provide excellent color accuracy to distinguish between tissues
- The ergonomic handle is off-center to maximize the light field for broader examinations

Ideal for use in:

ICU Rooms • ED / Trauma Departments • Procedure and Exam Rooms
Labor and Delivery Units • Cosmetic Surgery Center • Veterinary Surgery Center



80,000
LUX



4500K
CCT



95+
CRI



Effortless Positioning

Perform minor procedures and examinations with confidence.

Ergonomics


- The ultra slim design reduces clutter in small workplaces
- Lightweight long reach arms make positioning effortless
- 360 degrees of rotation for precise beam placement
- Easy access to lighting functions

Infection Control

- The removable handle can be sterilized to maximize your infection control efforts
- Fully sealed and gasketed housing provides a seamless finish for easy wipe down

Functionality

- Five levels of dimming allow you to adjust the light intensity with ease
- Endo mode provides a soft-ambient glow of blue light to the room to make it easier for clinicians to see the monitors and their surroundings during endoscopy procedures

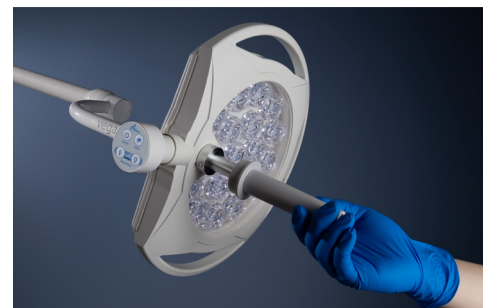

360 Degree
Rotation




Endo
Mode




50% - 100%
Dimming

The Vega 80 LED Series

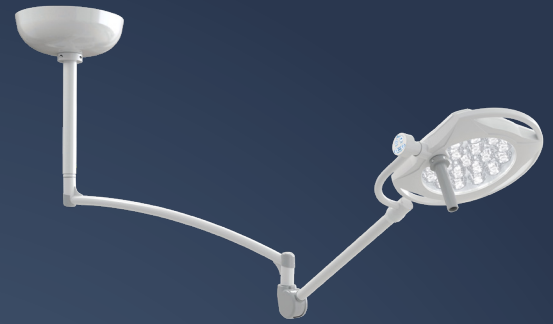

Ceiling Mount

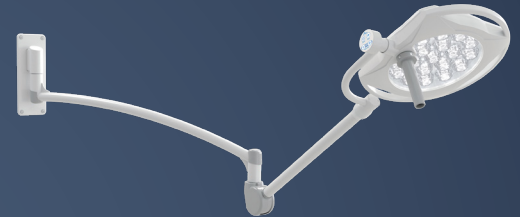

Wall Mount



Floor Stand

Lighthead Diameter	14" (35.6 cm)
Light Intensity	80,000 LUX @ 1m
Brightness Control	50% - 100%
Color Rendering Index	95+
R9	98
Lighthead Weight	7 lbs (3.3 kg)
Color Temperature (Standard)	4500K
Number of LEDs	18 LEDs
Light Field Diameter	7.5" (19 cm)
Depth of Illumination	94.5" (240 cm)
Electrical	100-240 VAC, 50/60 Hz, 0.5 amps
Power Supply	Direct Power Supply
LED Service Life	50,000 hours
Certifications	Conforms to: AAMI STD ES60601-1, IEC STD 60601-2-41; Certified to CSA C22.2 No. 60601-1
Warranty	Limited five-year warranty

Amico Lights Corporation
 122-B East Beaver Creek Road, Richmond Hill, ON L4B 1G6, Canada
 Toll Free Phone: 1.877.462.6426 | Toll Free Fax: 1.866.440.4986
 Tel: 905.764.0800 | Fax: 905.764.0862
 Email: info@amico.com | www.amico.com

Visit **amico.com** for more information
 about our **Vega 80 LED Exam Light**

