


GoLift Ultra



Amico

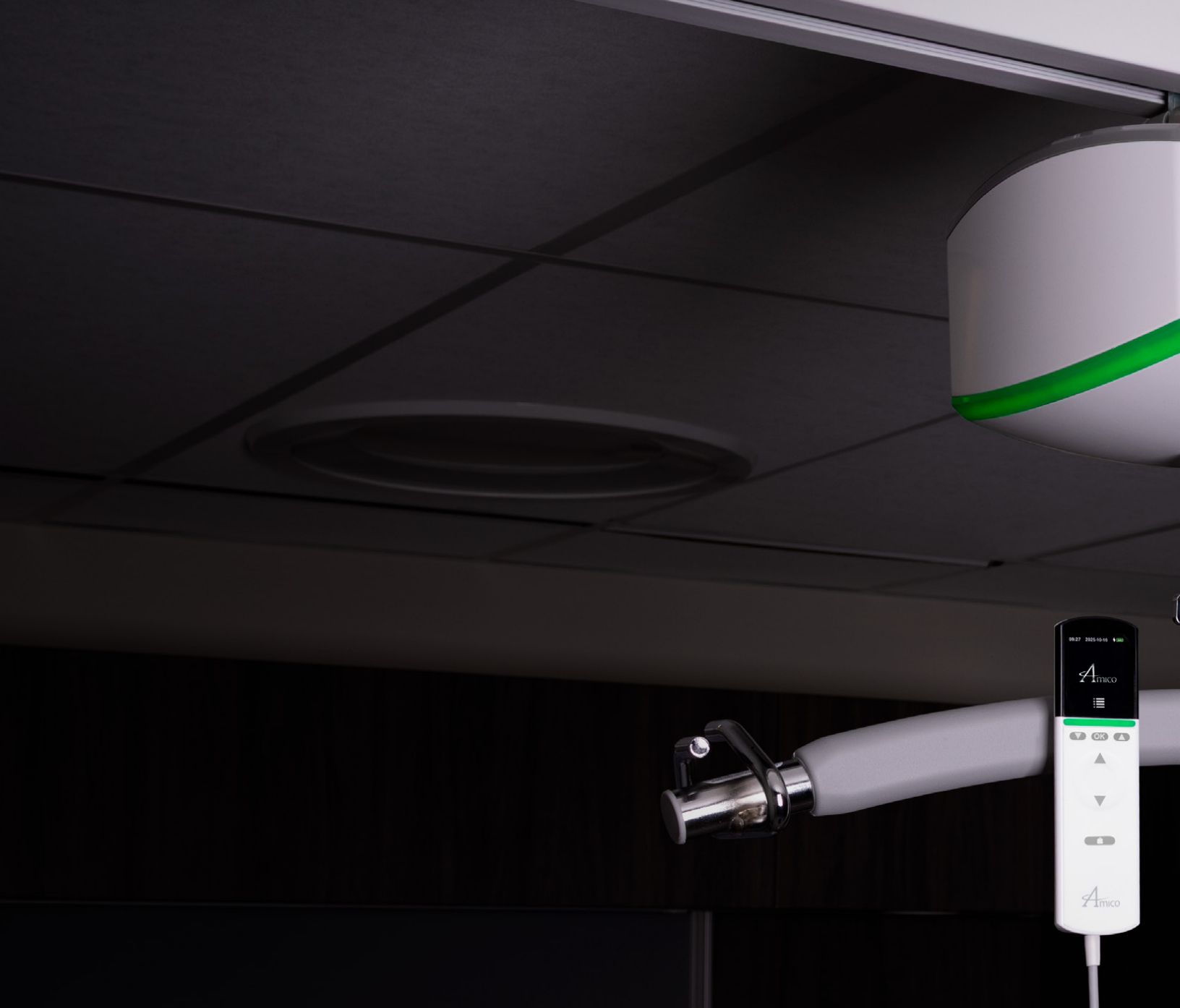


Amico Group of Companies

Amico Corporation, founded in 1974, designs, manufactures and markets a broad range of products for virtually every department in the health care facility. Our mission is to provide state-of-the-art products that facilitate the caregiver's role in providing the best outcomes for patients.

These products are manufactured in eight facilities in the U.S. and Canada, which occupy over **700,000** square feet.



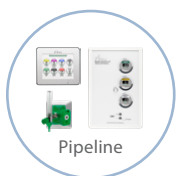


Pipeline



Architectural



OR

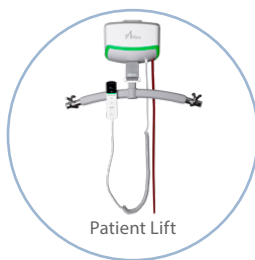

Patient Lift



ICU Booms




Lights



Patient Care



Beds

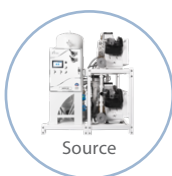

Source



Accessories



Diagnostic



SOT



Industrial


Amico Group of Companies - One Point of Contact. One Total Solution.

GoLift Ultra Features

- LED Indicator
- Integrated Scale with Memory Function
- Wired and Magnetic Digital Hand Control
- Hand Control with IP67
- User Friendly Interface (Battery Level, Warnings, Diagnostic Data and more)
- In-Rail Continuous or End-Stop Charging
- Wired Data Retrieval (Bluetooth or Wi-Fi)
- Quick Connect Carry Bar and Trolley
- Lightweight and Portable



Safe Patient Handling

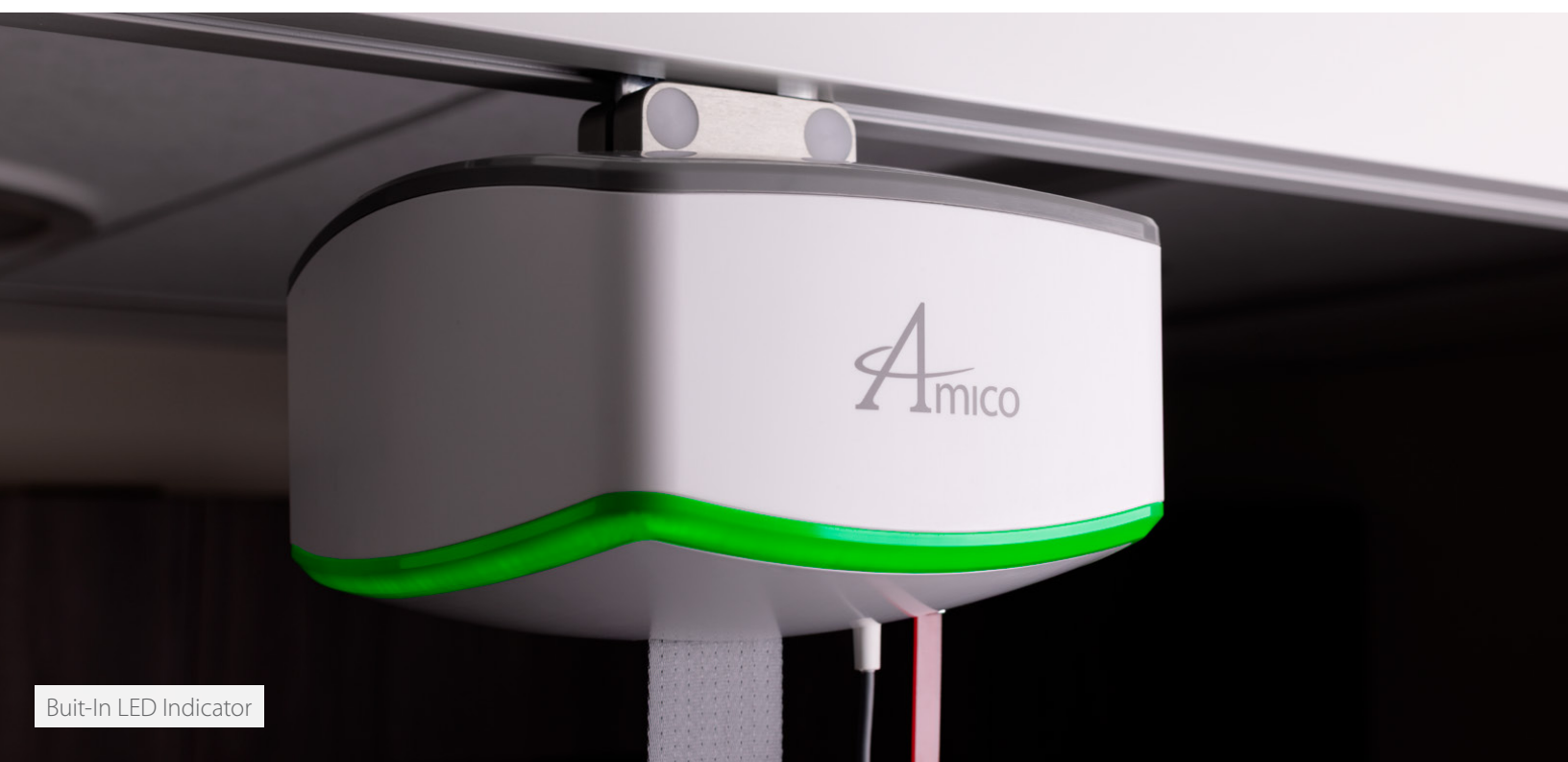

Built-In LED Indicator

GoLift Ultra

Safety Features

- Emergency Stop and Lowering Device
- Upper and Lower Limit Detection
- Slack Tape Sensor
- Free Fall Brake (Overspeed Governor)
- Low Battery and Dead Battery Alarms
- Soft Start and Stop
- Overload Protection
- Emergency Manual Lowering



Smallest Lift in the Market



Wired Magnetic Hand Control

GoLift Ultra Patient Lift Pendant (PLP)



Patient Lift Integration



Adapt Room to Patient



Side-to-Side Transfers and Reposition in Bed

GoLift Ultra Specifications



Maximum Lifting Speed

No load:	2.0 inches/second
GoLift Ultra450:	1.2 inch/second
GoLift Ultra700:	1.2 inch/second

Unit Weight

GoLift Ultra450:	14 lbs (6.3 kg)
GoLift Ultra700:	14 lbs (6.3 kg)

Safe Working Load (SWL)

GoLift Ultra450:	450 lbs (204 kg)
GoLift Ultra700:	700 lbs (318 kg)

** The maximum safe working load can be adjusted to a lower limit based on customer request*

Dimensions

GoLift Ultra450 / GoLift Ultra700:	
• Length:	8.9" (226.06 mm)
• Width:	8.9" (226.06 mm)
• Height:	6.3" (160.02 mm)
• Strap Length:	84" (2133.6 mm)

Number of Lifts Per Charge (Duty: 10/90)

25% of Strap at Midrange:

- 375 with 185 lbs (84 kg)
- 170 with 450 lbs (204 kg)
- 102 with 700 lbs (318 kg)

Charging time

2 - 4 hours

Lift Case

Flame Retardant ABS

Batteries

High Capacity, Nickel Metal Hydride (NiMH)

Standard: 2x 14.4V (3.3Ah)

With Over-Temperature and Over-Charge Protection

Battery Charger

Input:	100-240 VAC
Output:	36 VDC, 1.0A, 40W

Digital Hand Control

Hand Control with LCD Display (240x282 pixels)

Maximum Sound Level

<55 dB

Expected Service Life

10 years or 22,500 cycles

Approvals

Certified to:	ANSI/AAMI ES 60601-1:2012
	CSA-C22.2 No. 60601-1:2014
	UL94
Tested to:	ISO 10535:2021

Protection Class

IP67

Weigh Scale

Integrated Scale

Safety Features

Emergency Stop & Lowering Device
Upper and Lower Limit Detection
Slack Tape Sensor
Free Fall Brake (Overspeed Governor)
Low Battery and Dead Battery Alarms
Soft Start and Soft Stop
Overload Protection
Emergency Manual Lowering

Diagnostic/Usage Information

Number of Lifts:

- Total Number of Lifts
- Number of Heavy Lifts

Charging Statistics:


- Total Number of Charge Cycles

User Weight (If equipped)

Recent Faults

GoLift Ultra Patient Lift Pendant (PLP)

in the OR

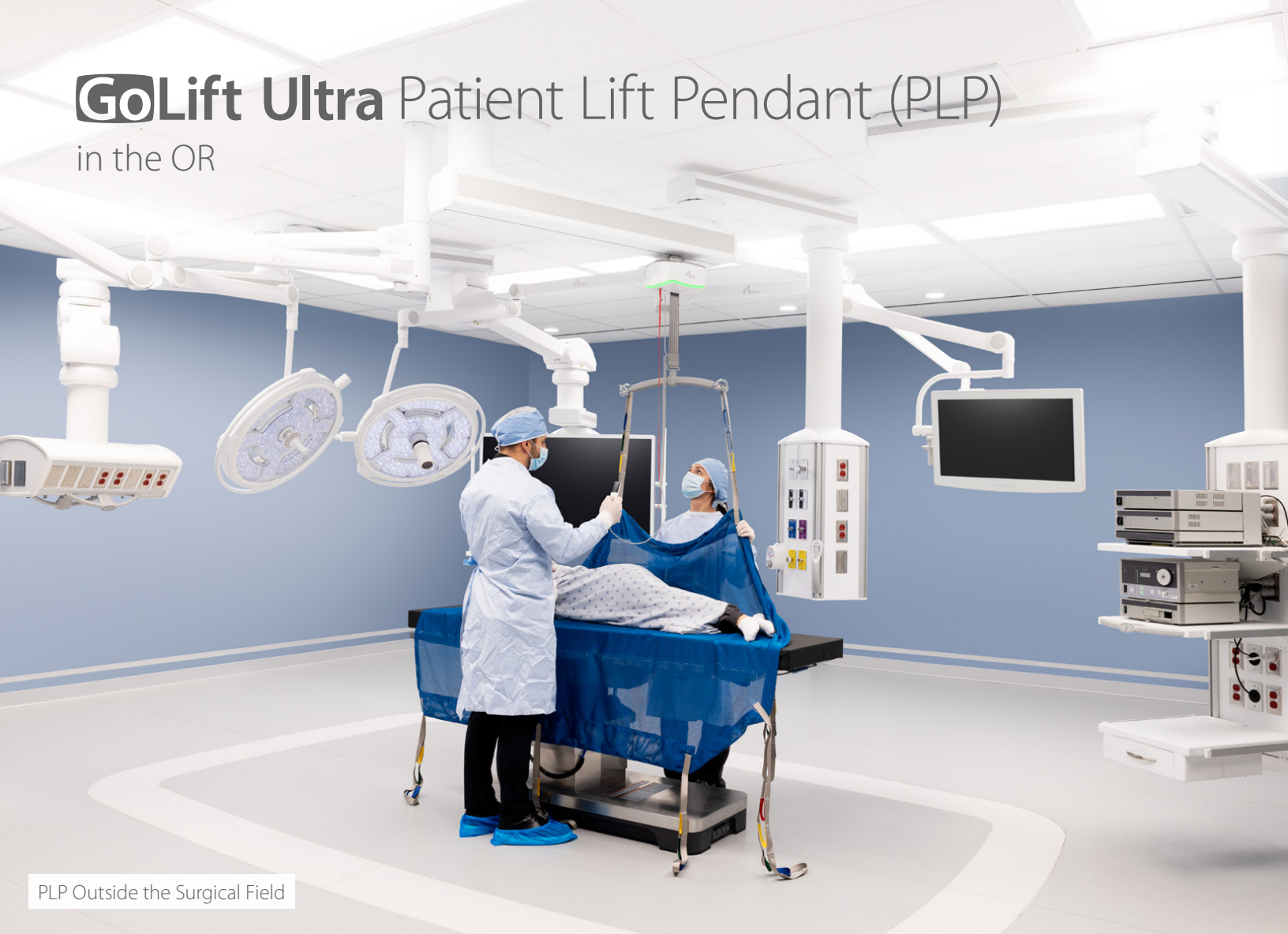

PLP Outside the Surgical Field

GoLift Ultra Patient Lift Pendant (PLP)

in the ICU



Reduced Risk of Injury to Patients and Caregivers


GoLift Ultra Headwall Integration



Hidden GoLift and Track



Easily Accessible

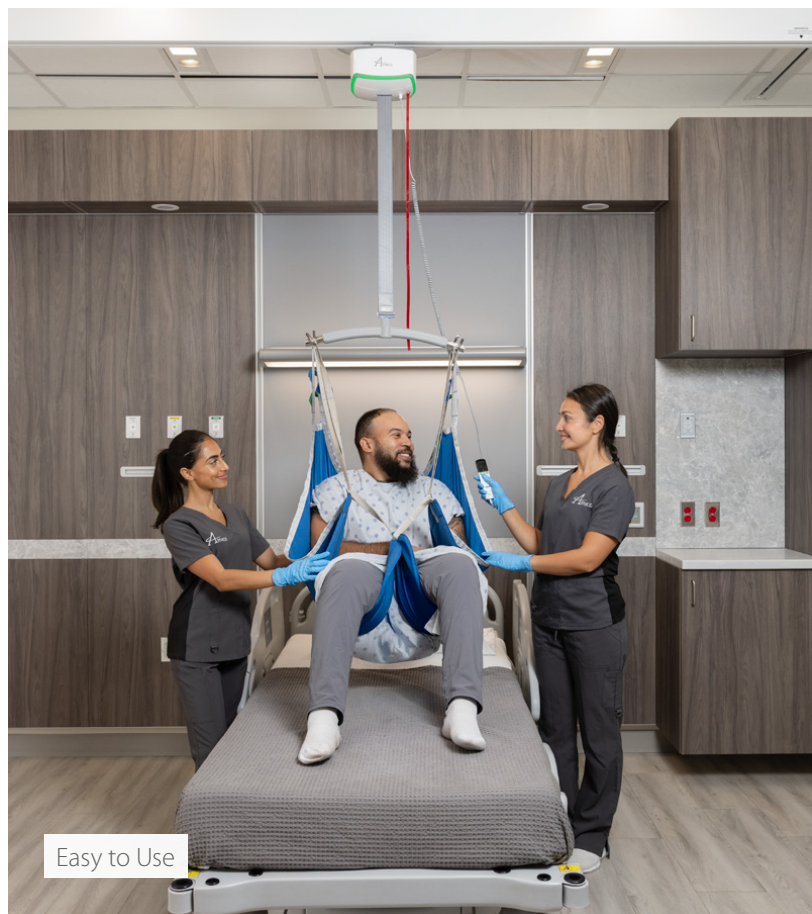

Easy to Use

GoLift Ultra Slings

We offer a wide selection of reusable and disposable slings. Refer to our Sling Brochure for more details.

Guidelines for Choosing your Sling Size

		S	M	L	XL
GoBasic		With Head Support Height: 29" (73.6 cm) Shoulder: 20" (50.8 cm) Hip: 16" (40.6 cm)	Height: 30" (76.2 cm) Shoulder: 24" (60.9 cm) Hip: 20" (50.8 cm)	Height: 31" (78.7 cm) Shoulder: 28" (71.1 cm) Hip: 24" (60.9 cm)	Height: 32" (81.2 cm) Shoulder: 32" (81.2 cm) Hip: 28" (71.1 cm)
	No Head Support	Height: 21" (53.3 cm) Shoulder: 20" (50.8 cm) Hip: 16" (40.6 cm)	Height: 22" (55.8 cm) Shoulder: 24" (60.9 cm) Hip: 20" (50.8 cm)	Height: 23" (58.4 cm) Shoulder: 28" (71.1 cm) Hip: 24" (60.9 cm)	Height: 24" (61 cm) Shoulder: 32" (81.2 cm) Hip: 28" (71.1 cm)
GoLong		42" x 78" (106.6 cm x 198.1 cm)		56" x 78" (142.2 cm x 198.1 cm)	
GoHygiene		Chest: 32" (81.2 cm)	Chest: 36" (91.4 cm)	Chest: 40" (101.6 cm)	Chest: 44" (111.7 cm)
GoComfort		Height: 23" (58.4 cm) Shoulder: 24" (60.9 cm) Hip: 20" (50.8 cm)	Height: 24" (60.9 cm) Shoulder: 26" (66 cm) Hip: 22" (55.8 cm)	Height: 25" (63.5 cm) Shoulder: 28" (71.1 cm) Hip: 24" (60.9 cm)	Height: 26" (66 cm) Shoulder: 30" (76.2 cm) Hip: 26" (66 cm)
GoActive		Chest: 40" (101.6 cm)	Chest: 44" (111.7 cm)	Chest: 48" (121.9 cm)	Chest: 52" (132 cm)
GoReposition		20" x 75" (51 cm x 190.5 cm)			
GoLimb		15" x 7" (38.1 cm x 17.8 cm)			
GoTurn		Width: 22" (55.8 cm)			
GoAmputee		Height: 28" (71.1 cm) Hip: 15" (38.1 cm)	Height: 33" (83.8 cm) Hip: 18" (45.7 cm)	Height: 38" (96.5 cm) Hip: 22" (55.9 cm)	Height: 43" (109.2 cm) Hip: 26" (66 cm)
GoWalk		Hip: 15" (38.1 cm)	Hip: 18" (45.7 cm)	Hip: 22" (55.9 cm)	Hip: 26" (66 cm)



GoTrack Systems

Configurations


- XY Gantry (H-Track)
- Straight Track
- Patient Lift Pendant
- Turntable (Manual or Power)
- Free-Standing Track
- Wall-Post Track
- Power Traverse Gantry

Features

- Light weight and high lifting
- Quick and easy
- Hidden end stops for safety
- Continuous charge or end stop charge

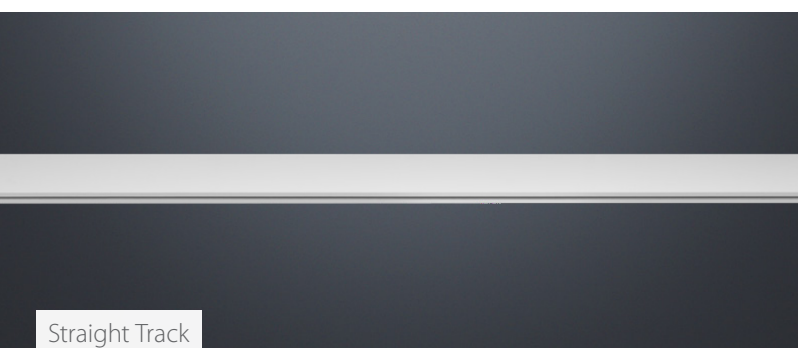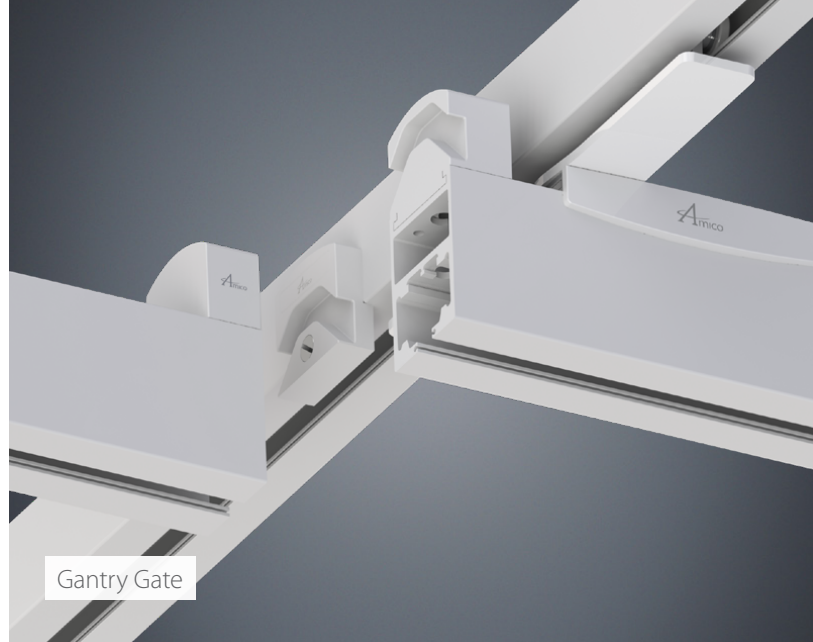

Anatomy of the **GoLift Ultra**

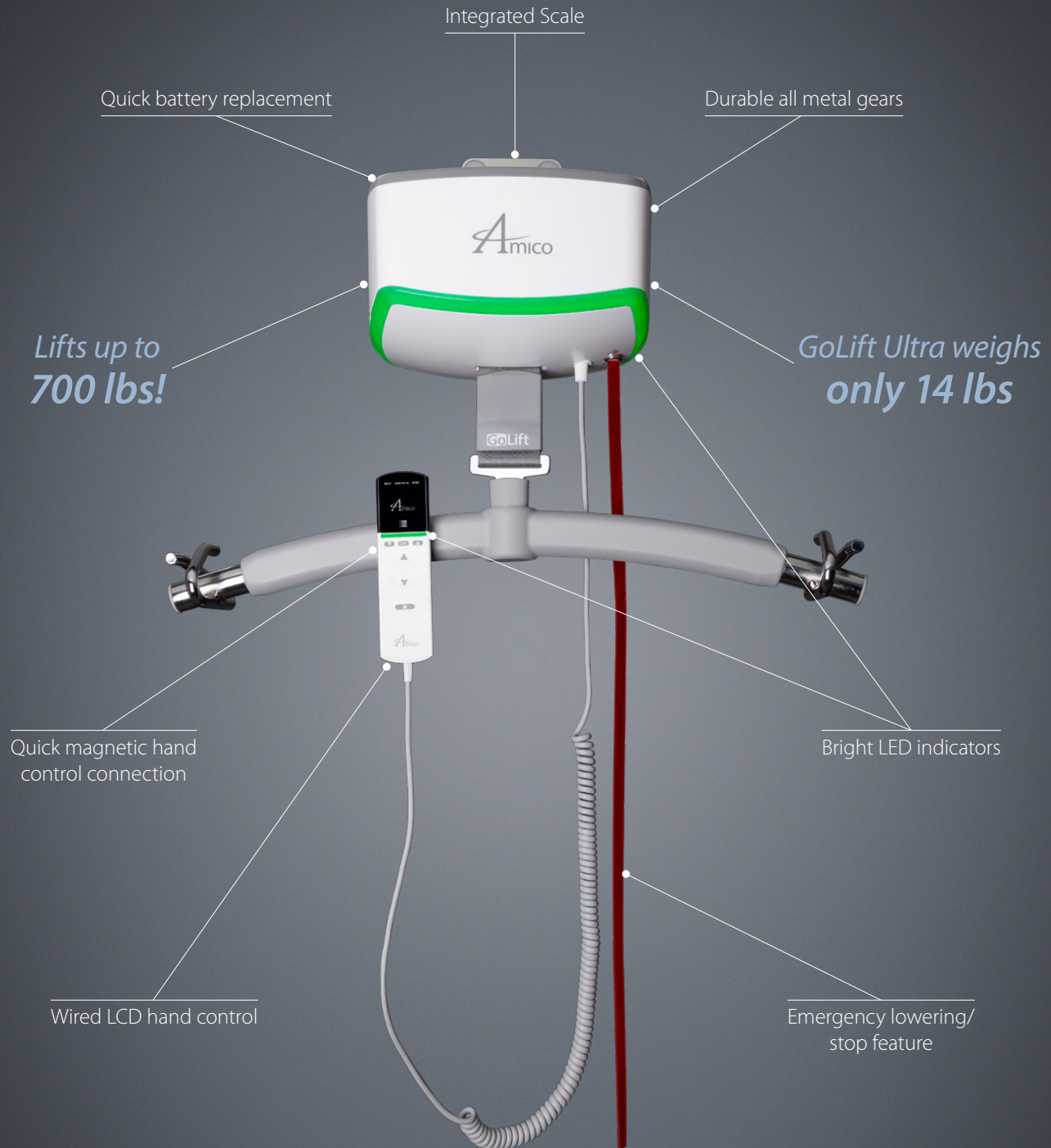

Amico Mobility Solutions Corporation | 122-B East Beaver Creek Road, Richmond Hill, ON L4B 1G6, Canada
600 Prime Place, Hauppauge, NY 11788, USA
Toll Free Tel: 1.877.462.6426 | Tel: 905.764.0800 | Fax: 905.764.0862
Email: info@amico.com | www.amico.com